

# ThevoTherapy - The therapy chair with mobility seat.

## Thevo - The mobility seat for children: Stability AND movement at the same time

The mobility seat now offers optimum support to the children and at the same time the possibility for small movements of the trunk. This finally creates the best possible conditions for the children to carry out well directed activities with their arms and hands.

### ThevoTherapy

- ✓ Individual adjustments by means of the flexible elements in the back.
- ✓ Firm seat as basis for an optimum seating posture.
- ✓ Height adjustable for an optimum use in everyday life.
- ✓ Growing with the child and easily adjustable to meet the individual needs.

### Therapy-frame

Comfortable indoor base, easy adjustment of height and seat tilt-in-space by gas spring (foot engaged). 4 swivelling castors, 2 of them lockable (rear).



#### Headrest

Adjustable in angle and height. The lateral neck support is individually adjustable.

#### Suspension unit

Adapts to the body contours and offers movement possibilities to the child.

#### Arm rests

Ergonomically-shaped arm rest, adjustable in width, height, and angle.

#### Gas spring

Back angle easily adjustable by gas spring.


#### Wattle interlaced yarns

Easy assembly of the accessories. No more need for slitting the fabric. Optimised for the re-use.

#### Anatomically shaped seat upholstery

Seat bone tuberosity padding with visco-elastic foam, low profile pommel.


### Accessories (selection, please visit our website to see the whole range)

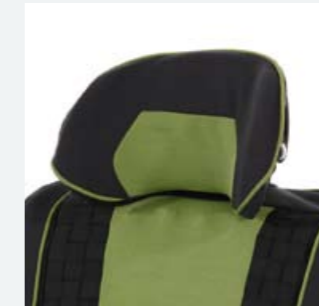


Headrest, shell: provides lots of lateral support, individually adaptable




Head rest, neck support: provides for stability by supporting the lower neck, individually adaptable

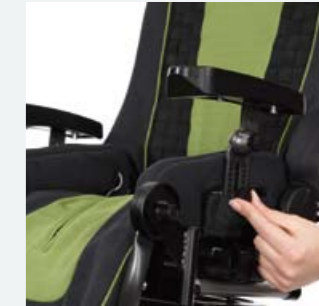


Arm rest

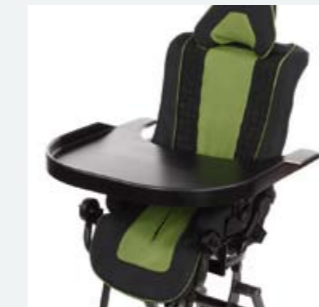

Tray

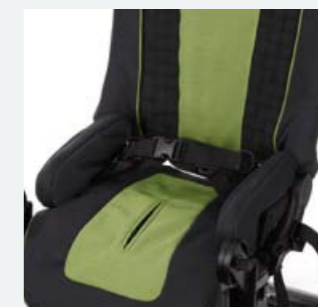


4-point pelvic belt with front pull

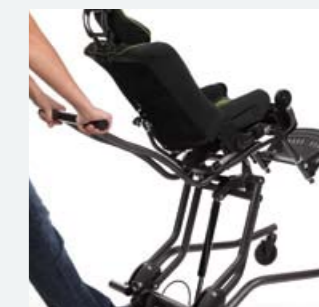

Push handle



Belt bow

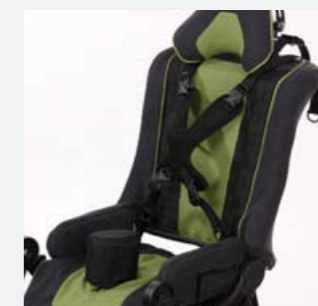

Butterfly/chest harness and abduction block

@ Further product info on the internet: [www.thomashilfen.com](http://www.thomashilfen.com)



## Technical data

ThevoTherapy	size 1	size 2
Seat depth	19 - 30 cm / 7.5 - 11.8"	27 - 38 cm / 10.6 - 15"
Seat width	19 - 31 cm / 7.5 - 12.2"	25 - 36 cm / 9.8 - 14.2"
Back height	47 - 68 cm / 18.5 - 26.8"	64 - 78 cm / 25.2 - 30.7"
Lower leg length	17 - 30 cm / 6.7 - 11.8"	25 - 36 cm / 9.8 - 14.2"
Knee angle (adjustable in 15° steps)	90° to 180°	90° to 180°
Hip angle	90° to 122°	90° to 122°
Foot angle	78° to 100°	78° to 100°
Size of foot rest (w x l)	28.5 x 18 cm / 11.2 x 7"	33 x 21 cm / 12.9 x 8.3"
Overall dimensions seat (w x h x l)	38 x 78 x 54 cm / 15 x 30.7 x 21.3"	45 x 101 x 66 cm / 17.7 x 39.8 x 26"
Folding dimensions seat (w x h x l)	38 x 41 x 57 cm / 15 x 16.1 x 22.4"	45 x 48 x 74 cm / 17.7 x 18.9 x 29.1"
Carrying capacity seat	35 kg / 77.2 lb	40 kg / 88.2 lb
Seat tilt-in-space	-4° to +31°	-4° to +31°
Seat height with seat unit	26 - 66 cm / 10.2 - 26"	26 - 66 cm / 10.2 - 26"
Overall dimensions chassis (w x h x l)	54 x 59 x 73 cm / 21.3 x 23.2 x 28.7"	54 x 59 x 73 cm / 21.3 x 23.2 x 28.7"
Folding dimensions chassis (w x h x l)	54 x 24 x 73 cm / 21.3 x 9.5 x 28.7"	54 x 24 x 73 cm / 21.3 x 9.5 x 28.7"
Push bar height (accessory)	80 - 102 cm / 31.5 - 40.2"	80 - 102 cm / 31.5 - 40.2"
Adapter height	20 - 59 cm / 7.9 - 23.2"	20 - 59 cm / 7.9 - 23.2"
Size of swivel wheels (front)	100 mm / 3.9"	100 mm / 3.9"
Wheel size (rear with brake)	100 mm / 3.9"	100 mm / 3.9"
Carrying capacity chassis	50 kg / 110.2 lb	50 kg / 110.2 lb
Weight	19.9 kg / 43.9 lb	21.4 kg / 47.2 lb



ThevoTherapy

Thevo  
by Thomashilfen



Nursing scientific consultation by



Checked by



Item code: 95732, status: 10/2009, subject to technical alterations

Thevo  
by Thomashilfen

Walkmühlenstraße 1  
D - 27432 Bremervörde · Germany  
Phone: +49 (0)47 61 / 886-63 or -68  
Fax: +49 (0)47 61 / 886-19  
www.thevo-division.com



# Thevo - The mobility seat