



**NEW
size 2**

Swiftly

The lightweight rehab strollers
for both small & big kids



Simple use, light & swift: Here comes **Swiftly** ...

“ A stable, comfortable rehab stroller with a fixed seat; it can withstand the many stresses of day-to-day life and is quickly ready for use. Now available in two sizes for growing children. ”

Comfortable sitting for the child: the washable seat cover is soft and breathable – preventing perspiration and providing a thoroughly cozy feeling



Back angle can be adjusted for positioning and relaxed seating



Fold lock (only for size 2) keeps everything securely folded together.



size 2



Small folding size (even size 2) – Especially quick and easy to fold with a flick of the wrist.



size 1

With the spacious NEW Swiftly size 2

- big kids that weigh up to 50 kg / 110.2 lb
 - and are up to 150 cm / 59.1" tall
- can be cared for with a rehab stroller

Seat, back and footrest provide safe and firm support – a convenient place to sit for the child

H-belt with padded harness as a standard feature



Seat minimizer for smaller children (accessory) – grows with your child for a very long time



Seat depth adjustable without tools – even while the child is sitting in the rehab stroller (only for size 1)



Foldable and height-adjustable footrest – Children can get in and out of the stroller independently.

Transit tie-down-kit for safe transport in a vehicle (available for both sizes, standard for size 2)



Swiftly & Swiftly 2

The lightweight rehab strollers for both small & big kids

Swiftly – Rehab stroller



All details online:
www.thomashilfen.us/swifty1



or .../swifty2



Accessories (selection)



Sleeping bag,
 padded (summer)
 or woven fur (winter)
(shown on size 1)



Sun and rain cover
 with leg blanket
(shown on size 1)



Umbrella
(shown on size 1)




Basket
(shown on size 1)



Butterfly chest harness
 (flexible), abduction block
 and grip rail
(shown on size 2)



Head pillow, lateral trunk
 support with chest belt and
 2-point pelvic harness
(shown on size 2)

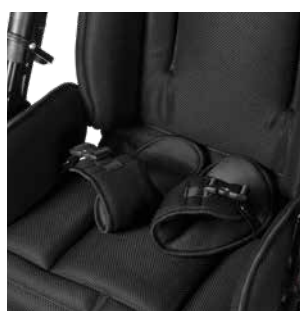

Soft pelvic harness
(shown on size 2)



Tip assist
(shown on size 2)

You will find the total overview of all rehab accessories
 on our website: www.thomashilfen.us

Technical data	Swiftly 1	Swiftly 2
Seat depth (with seat minimizer)	22 (17 cm) - 28.5 cm / 8.7 (6.9) - 11.2"	37.5 cm (28 cm) / 14.8" (11")
Seat width (with seat minimizer)	34 cm (23 cm) / 13.4" (9.1")	39.5 cm (29.5 cm) / 15.6" (11.6")
Back height	62 cm / 24.4"	70 cm / 27.6"
min. shoulder height with seat minimizer	26 cm / 10.2"	30 cm / 11.8"
Lower leg length (with seat minimizer)	16 (12 cm) - 33 cm / 6.3 (4.7") - 13"	21 (21 cm) - 37 cm / 8.3 (8.3 ") - 14.6"
Hip angle	90° to 115°	90° to 100°
Foot angle	90°	90°
Seat tilt-in-space fix	+15°	+20°
Size of foot rest (w x d)	31 x 18 cm / 12.2 x 7.1"	31 x 20 cm / 12.2 x 7.9"

Technical data	Swiftly 1	Swiftly 2
Push bar height	80 - 119 cm / 31.5 - 46.9"	82 - 124 cm / 32.3 - 48.8"
Wheel size swivel wheels (front)	18 cm / 7.1"	18 cm / 7.1"
Wheel size (rear)	25 cm / 9.8"	25 cm / 9.8"
Overall dimensions (l x w x h)	61 x 102 x 98 cm / 24 x 40.2 x 38.6"	67 x 113 x 98 cm / 26.4 x 44.5 x 38.6"
Folded dimensions (l x w x h)	61 x 38 x 73 cm / 24 x 15 x 28.7"	67 x 49 x 79 cm / 26.4 x 19.3 x 31.1"
Weight	12.4 kg / 27.3 lb	18.2 kg / 40.1 lb
max. user weight	35 kg / 77.2 lb	50 kg / 110.2 lb



Thomashilfen North America (ExoMotion, LLC) · 7936 Occidental Ave S · Seattle WA 98108 · USA
 Phone: 866.870.2122 (toll free) · Fax: 866.870.0801 (toll free) · info@thomashilfen.com · www.thomashilfen.us